

Antibody Data

Product SKU:	AGEL1367	Clone:	AFS98
Applications:	FCM		
Reactivity:	Mouse		

Important Note:


Centrifuge before opening to ensure complete recovery of vial contents.

Product Information:

Alternate Names: Macrophage colony-stimulating factor 1 receptor;CSF-1 receptor (EC:2.7.10.1);CSF-1-R;CSF-1R;M-CSF-R;Proto-oncogene c-Fms;CD115;Csf1r;Csfmr;Fms;
Uniprot ID: P09581

Background: CSF-1R, also known as CD115 and M-CSFR, is a single-pass type I membrane protein and member of the platelet-derived growth factor receptor family. This c-fms (Fms proto-oncogene) gene product's natural ligands include M-CSF and IL-34. Structural studies of CD115 have described an Ig-like extracellular domain, a transmembrane domain, an intracellular juxtamembrane domain, a split tyrosine kinase domain, and a C-terminal tail receptor. Receptor activation induces homodimerization in addition to phosphorylation and ubiquitination of intracellular residues. CD115 directly influences tissue macrophage and osteoclast differentiation and proliferation. It is expressed on monocytes/macrophages, peritoneal exudate cells, plasmacytoid and conventional dendritic cells, and osteoclasts.

Form:	Liquid
Conjugation:	APC
Size:	50 Tests, 100 Tests, 200 Tests
Host Species:	Rat
Isotype:	Rat IgG2a, κ

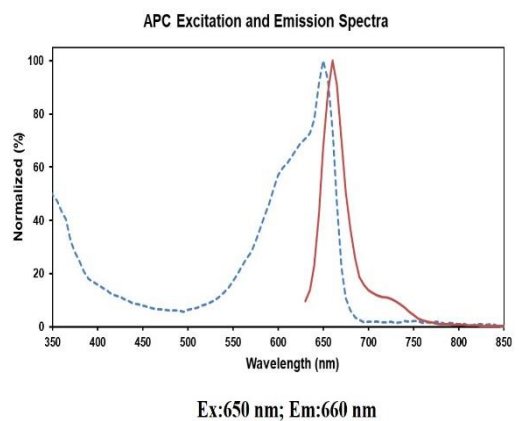

Isotype Control:	APC Rat IgG2a, κ Isotype Control[2A3] [Product AGEL1367]
Storage Buffer:	Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.
Shipping:	Biological ice pack at 4°C

Stability & Storage: Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial contents. This product is guaranteed up to one year from purchase.

Recommended Usage: Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.
