

**Antibody Data**

<b>Product SKU:</b>	<b>AGEL1630</b>	<b>Clone:</b>	<b>29F.1A12</b>
<b>Applications:</b>	<b>FCM</b>		
<b>Reactivity:</b>	<b>Mouse</b>		

**Important Note:**

Centrifuge before opening to ensure complete recovery of vial contents.

**Product Information:**

**Alternate Names:** PD-1; Programmed Death-1;

**Uniprot ID:** Q02242

**Background:** CD279, also known as programmed death-1 (PD-1), is a 50-55 kD glycoprotein belonging to the CD28 family of the Ig superfamily. PD-1 is expressed on activated splenic T and B cells and thymocytes. It is induced on activated myeloid cells as well. PD-1 is involved in lymphocyte clonal selection and peripheral tolerance through binding its ligands, B7-H1 (PD-L1) and B7-DC (PD-L2). It has been reported that PD-1 and PD-L1 interactions are critical to positive selection and play a role in shaping the T cell repertoire. PD-L1 negative costimulation is essential for prolonged survival of intratesticular islet allografts.


**Form:** Liquid

**Conjugation:** PE/Cyanine 5.5

**Size:** 50 Tests, 100 Tests, 200 Tests

**Host Species:** Rat

**Isotype:** Rat IgG2a, κ

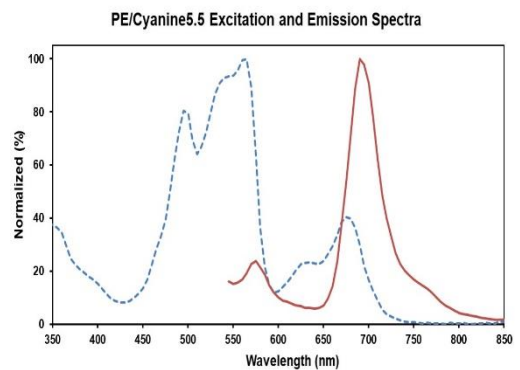

Ex:495;565;675 nm; Em:690 nm

**Isotype Control:** PE/Cyanine5.5 Rat IgG2a, κ Isotype Control[2A3] [Product AGEL1630]

**Storage Buffer:** Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

**Shipping:** Biological ice pack at 4°C

**Stability & Storage:** Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial contents. This product is guaranteed up to one year from purchase.

**Recommended Usage:** Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu$ L of antibody per test (million cells in 100  $\mu$ L staining volume or per 100  $\mu$ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.