

#### **Product Datasheet**

# PE/Cyanine7 Anti-Human CD45RO Antibody [UCHL1]

Catalogue Code: AGEL1701

### **Antibody Data**

Product SKU: AGEL1701 Clone: UCHL1

Applications: FCM

Reactivity: Human

## **Important Note:**

Centrifuge before opening to ensure complete recovery of vial contents.

#### **Product Information:**

Alternate Names: Uniprot ID: -

**Background**: CD45RO is a 180 kD single chain membrane glycoprotein. It is a splice variant of tyrosine

phosphatase CD45, lacking the A, B, and C determinants. The CD45RO isoform is expressed on activated and memory T cells, some B cell subsets, activated monocytes/macrophages, and granulocytes. CD45RO enhances both T cell receptor and B cell receptor signaling mediated activation. CD45 and its isoforms non-covalently associate with lymphocyte phosphatase-associated phosphoprotein (LPAP) on T and B lymphocytes. CD45 has been reported to be associated with several other cell surface antigens including CD1, CD2, CD3, and CD4. CD45 has also been reported to bind galectin-1 and CD22. CD45 isoform expression can change in response to cytokines.


Form: Liquid

**Conjugation:** PE/Cyanine 7

Size: 20 Tests, 100 Tests, 200 Tests

Host Species: Mouse

**Isotype:** Mouse IgG2a, κ

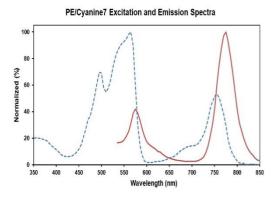

Ex:495;565;755 nm; Em:775 nm

**Isotype Control:** PE/Cyanine7 Mouse IgG2a, κ Isotype Control[C1.18.4] [Product AGEL1701]

**Storage Buffer:** Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

**Shipping:** Biological ice pack at 4°C



Stability & Storage:

Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial contents. This product is guaranteed up to one year from purchase.

Recommended Usage:

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu$ L of antibody per test (million cells in 100  $\mu$ L staining volume or per 100  $\mu$ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.