

**GenieFluor 488 Anti-Human CD45RO  
Antibody [UCHL1]  
Catalogue Code: AGEL1704**

**Antibody Data**

<b>Product SKU:</b>	<b>AGEL1704</b>	<b>Clone:</b>	<b>UCHL1</b>
<b>Applications:</b>	<b>FCM</b>		
<b>Reactivity:</b>	<b>Human</b>		

**Important Note:**

Centrifuge before opening to ensure complete recovery of vial contents.

**Product Information:**

**Alternate Names:** -

**Uniprot ID:** -

**Background:** CD45RO is a 180 kD single chain membrane glycoprotein. It is a splice variant of tyrosine phosphatase CD45, lacking the A, B, and C determinants. The CD45RO isoform is expressed on activated and memory T cells, some B cell subsets, activated monocytes/macrophages, and granulocytes. CD45RO enhances both T cell receptor and B cell receptor signaling mediated activation. CD45 and its isoforms non-covalently associate with lymphocyte phosphatase-associated phosphoprotein (LPAP) on T and B lymphocytes. CD45 has been reported to be associated with several other cell surface antigens including CD1, CD2, CD3, and CD4. CD45 has also been reported to bind galectin-1 and CD22. CD45 isoform expression can change in response to cytokines.


**Form:** Liquid

**Conjugation:** Genie Fluor488

**Size:** 20 Tests, 100 Tests, 200 Tests

**Host Species:** Mouse

**Isotype:** Mouse IgG2a, κ

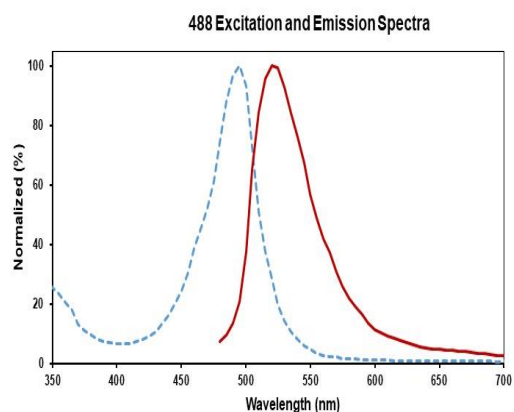

**Isotype Control:** Genie Fluor 488 Mouse IgG2a, κ Isotype Control[C1.18.4] [Product AGEL1704]

**Storage Buffer:** Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

**Shipping:** Biological ice pack at 4°C

**Stability & Storage:** Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial contents. This product is guaranteed up to one year from purchase.

**Recommended Usage:** Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu$ L of antibody per test (million cells in 100  $\mu$ L staining volume or per 100  $\mu$ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.