

#### **Product Datasheet**

# **GenieFluor 647 Anti-Human CD11b Antibody [ICRF44]**

Catalogue Code: AGEL1758

### **Antibody Data**

Product SKU: AGEL1758 Clone: ICRF44

Applications: FCM

Reactivity: Human

## **Important Note:**

Centrifuge before opening to ensure complete recovery of vial contents.

#### **Product Information:**

Alternate Names: Integrin alpha-M;ltgam;CD11 antigen-like family member B;CR-3 alpha chain;Leukocyte

adhesion receptor MO1;CD11b;

Uniprot ID: P11215

**Background**: CD11b is a 165-170 kD type I transmembrane glycoprotein also known as αMintegrin,

Mac-1, CR3, and C3biR.  $\dot{\text{CD}}$ 11b non-covalently associates with integrin  $\beta2(\text{CD18})$  and is expressed on granulocytes, monocytes/macrophages, dendritic cells, NK cells, and subsets of T and B cells. CD11b/CD18 is critical for the transendothelial migration of monocytes and neutrophils. It is also involved in granulocyte adhesion, phagocytosis, and neutrophil activation. CD11b/CD18 interacts with ICAM-1 (CD54), ICAM-2 (CD102),

ICAM-4, CD14, CD23, heparin, iC3b, fibrinogen, and factor X.


Form: Liquid

**Conjugation:** Genie Fluor647

Size: 20 Tests, 100 Tests, 200 Tests

Host Species: Mouse

**Isotype:** Mouse lgG1,  $\kappa$ 

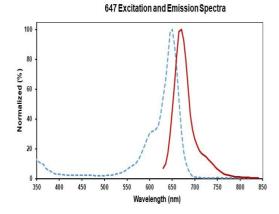

**Isotype Control:** Genie Fluor 647 Mouse IgG1, κ Isotype Control[MOPC-21] [Product AGEL1758]

Storage Buffer: Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

**Shipping:** Biological ice pack at 4°C



**Stability & Storage:** Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to

light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial

contents. This product is guaranteed up to one year from purchase.

Recommended Usage:

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu$ L of antibody per test (million cells in 100  $\mu$ L staining volume or per 100  $\mu$ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.