

**GenieFluor 488 Anti-Mouse CD117  
Antibody [2B8]**

Catalogue Code: AGEL2179

**Antibody Data**

<b>Product SKU:</b>	<b>AGEL2179</b>	<b>Clone:</b>	<b>2B8</b>
<b>Applications:</b>	<b>FCM</b>		
<b>Reactivity:</b>	<b>Mouse</b>		

**Important Note:**

Centrifuge before opening to ensure complete recovery of vial contents.

**Product Information:**

**Alternate Names:** Mast/stem cell growth factor receptor Kit;Kit;SCFR;Proto-oncogene c-Kit;Tyrosine-protein kinase Kit;CD117;

**Uniprot ID:** P05532

**Background:** CD117 is a 145 kD immunoglobulin superfamily member also known as c-Kit and stem cell factor receptor (SCFR). It is a transmembrane tyrosine-kinase receptor that binds the c-Kit ligand (also known as steel factor, stem cell factor, and mast cell growth factor). CD117 is expressed on hematopoietic stem cells (including multipotent hematopoietic stem cells, progenitors committed to myeloid and/or erythroid lineages, and T and B cell precursors), mast cells, and acute myeloid leukemia (AML) cells. CD117 interaction with its ligand is critical for the development of hematopoietic stem cells.


**Form:** Liquid

**Conjugation:** Genie Fluor488

**Size:** 50 Tests, 100 Tests, 200 Tests

**Host Species:** Rat

**Isotype:** Rat IgG2b, κ

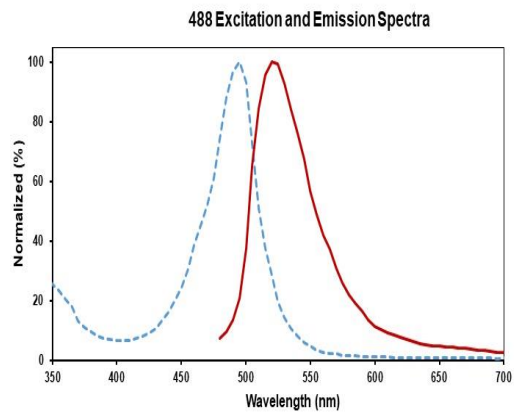

**Isotype Control:** Genie Fluor 488 Rat IgG2b, κ Isotype Control[LTF-2] [Product AGEL2179]

**Storage Buffer:** Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

**Shipping:** Biological ice pack at 4°C

**Stability & Storage:** Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial contents. This product is guaranteed up to one year from purchase.

**Recommended Usage:** Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu$ L of antibody per test (million cells in 100  $\mu$ L staining volume or per 100  $\mu$ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.