

#### **Product Datasheet**

# PerCP Anti-Human CD16 Antibody [3G8]

Catalogue Code: AGEL2480

## **Antibody Data**

Product SKU: AGEL2480 Clone: 3G8

Applications: FCM

Reactivity: Human

### **Important Note:**

Centrifuge before opening to ensure complete recovery of vial contents.

## **Product Information:**

**Alternate Names:** FcyRIII;Fc gamma receptor;Fc gamma receptor 3;

**Uniprot ID:** P08637 O75015

Background: CD16 is a 60 kD highly glycosylated protein. It is a member of the Ig superfamily and is

also known as B7-1, B7, and Ly-53. CD16 is constitutively expressed on dendritic cells and monocytes/macrophages, and inducibly expressed on activated B and T cells. The ligation of CD28 on T cells with CD16 and CD86 (B7-2) on antigen presenting cells (such as dendritic cells, macrophages, and B cells) elicits co-stimulation of T cells resulting in enhanced cell activation, proliferation, and cytokine production. CD16 appears to be expressed later in the immune response than CD86. CD16 can also bind to CD152, also

known as CTLA-4, to deliver an inhibitory signal to T cells.


Form: Liquid

Conjugation: PerCP

Size: 20 Tests, 100 Tests, 200 Tests

Host Species: Mouse

**Isotype:** Mouse IgG1, κ

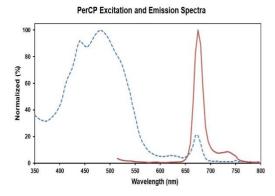

Ex:440;480;675 nm; Em:675 nm

**Isotype Control:** PerCP Mouse IgG1, κ Isotype Control[MOPC-21] [Product AGEL2480]

**Storage Buffer:** Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

**Shipping:** Biological ice pack at 4°C



**Stability & Storage:** Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to

light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial

contents. This product is guaranteed up to one year from purchase.

Recommended Usage:

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu$ L of antibody per test (million cells in 100  $\mu$ L staining volume or per 100  $\mu$ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.