

Product Datasheet

FITC Anti-Mouse CD200 Antibody [OX-

Catalogue Code: AGEL3385

Antibody Data

Product SKU: AGEL3385 Clone: OX-90

Applications: FCM

Reactivity: Mouse

Important Note:

Centrifuge before opening to ensure complete recovery of vial contents.

Product Information:

Alternate Names: OX-2 membrane glycoprotein; OX-2; MRC;

Uniprot ID: O54901

Background: CD200 (OX-2 antigen) is a type-1 membrane glycoprotein containing two extracellular Ig-

like domains. CD200 a highly conserved type I membrane glycoprotein that is expressed on a variety of cell types including thymocytes, some T cells, endothelial and follicular dendritc cells, B cells, and brain tissue (neurons); but not on NK cells, granulocytes, monocytes, or macrophages. CD200 costimulates T cell proliferation. It may regulate myeloid cell activity in a variety of tissues. CD200 is the ligand for CD200 receptor (CD200R). The CD200 Receptor is restricted to myeloid cells, and it is believed that its engagement with CD200 results in inhibition and/or downregulation of myeloid cell activity. Blocking of CD200/CD200R interactions decreases myeloid cell inhibitory thresholds

which results in enhanced immune activation.


Form: Liquid

Conjugation: **FITC**

Size: 50 Tests, 100 Tests, 200 Tests

Host Species: Rat

Isotype: Rat IgG2a, ĸ

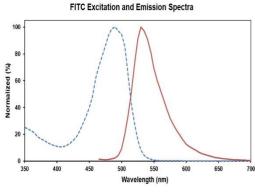

Ex:490 nm; Em:530 nm

FITC Rat IgG2a, κ Isotype Control[2A3] [Product AGEL3385] **Isotype Control:**

Storage Buffer: Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping: Biological ice pack at 4°C



Stability & Storage:

Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial contents. This product is guaranteed up to one year from purchase.

Recommended Usage:

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.